

Juni 2017

Oscillating cast cutter standard 500W, new improved heavy duty pattern

for all kind of Plaster of Paris bandages (P.O.P. bandages) and synthetic casts, even big and thick walled bandages
- cast dust aspiration available upon request



30.216.22 Oscillating cast cutter 500 W “**heavy duty**“, 230 V~, 50/60Hz and 500 Watt, with variable control of oscillations, **very silent motor** with an ergonomic housing, can be used for big and even thick walled bandages such as bandages for the thoracic and pelvic area.

The variation of oscillations per minute is from 3500 osc./min up to 24000 osc./min.

Ready for use complete with:

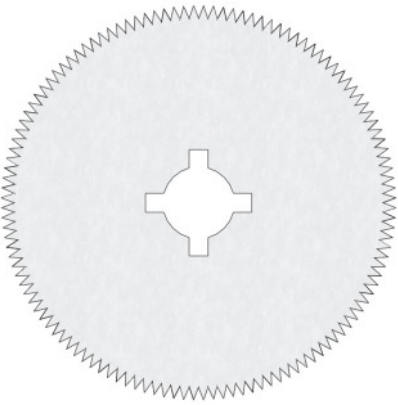
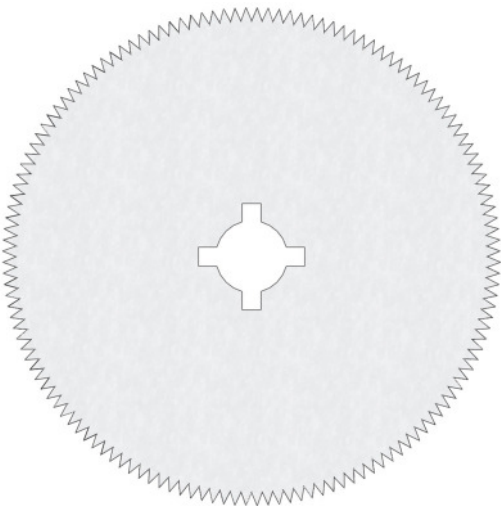
1 piece **30.211.50** saw blade 50 mm „**hard**“ for all **synthetic casts and Plaster of Paris bandages (P.O.P. bandages)**

1 piece **30.211.65** saw blade 65 mm „**hard**“ for all **synthetic casts and Plaster of Paris bandages (P.O.P. bandages)**

1 pair **30.212.02** wrench to exchange the saw blade

Technical details:	Code number 30.216.22
Current :	230 V~ 50/60 Hz
Length of cord:	approx. 5 m with EURO-plug
Power consumption:	500 Watt
Tested in accordance to:	EN60601-1 and -2, IEC 60601-1 and -2
Protection class:	II in accordance to VDE 0740 and type BF
Insulation:	double protection insulation
Switch for overloaded motor :	on the rear side of the machine
Support of motor:	plaster dust protected, with reinforced armature spindle
Weight:	1.500 gr. (with cord)

Saw blades for oscillating cast cutter standard 500 W

 <p style="text-align: center;">Illustration scale: 1: 1</p>	<p>30.217.50 Saw blade 50 mm „smooth cut” with special black coating for all synthetic bandages and Plaster of Paris (P.O.P) bandages</p> <p>30.211.50 Saw blade 50 mm “hard” for synthetic and P.O.P bandages</p> <p>30.210.50 Saw blade 50 mm for Plaster of Paris (P.O.P.) bandages</p>
 <p style="text-align: center;">Illustration scale: 1: 1</p>	<p>30.217.65 Saw blade 65 mm „smooth cut” with special black coating for all synthetic bandages and Plaster of Paris (P.O.P) bandages</p> <p>30.211.65 Saw blade 50 mm “hard” for synthetic and P.O.P bandages</p> <p>30.210.65 Saw blade 65 mm for Plaster of Paris (P.O.P.) bandages</p>